

Analysis Report

Official laboratory analysis summary for the submitted sample and associated quality-control review.

SAMPLE AOD-9604	RECEIVED DATE Jul 13, 2026	ANALYSIS DATE Jul 15, 2026	REPORT GENERATED Jul 17, 2026
STRENGTH	5mg	MANUFACTURER	PeptiCoreAminos
BATCH NUMBER	PC-AD5-0730U	LAB CODE	982-1
CLIENT	www.peptiCoreAminos.net		

SAMPLE INFORMATION

AOD-9604

5MG

FORM White lyophilized powder in a glass vial

SAMPLE SUBMISSION Sample provided by customer

BATCH PC-AD5-0730U

CAP / CRIMP COLOR Transparent / copper

RECEIVED DATE Jul 13, 2026




SAMPLE IMAGE

ANALYTICAL SUMMARY

IDENTITY	AOD-9604
PURITY	99.783%
QUANTITY	5.20mg
BATCH	PC-AD5-0730U
MANUFACTURER	PeptiCoreAminos

RP-HPLC-UV CHROMATOGRAM (220 NM)
Detection: UV 220 nm | Runtime: 12.0 min

Sample ID: 800-9684
Report ID: 578-2024-90019
Method: RP-HPLC-UV Method for Peptide Analysis
Detector: UV 220 nm | Runtime: 12.0 min

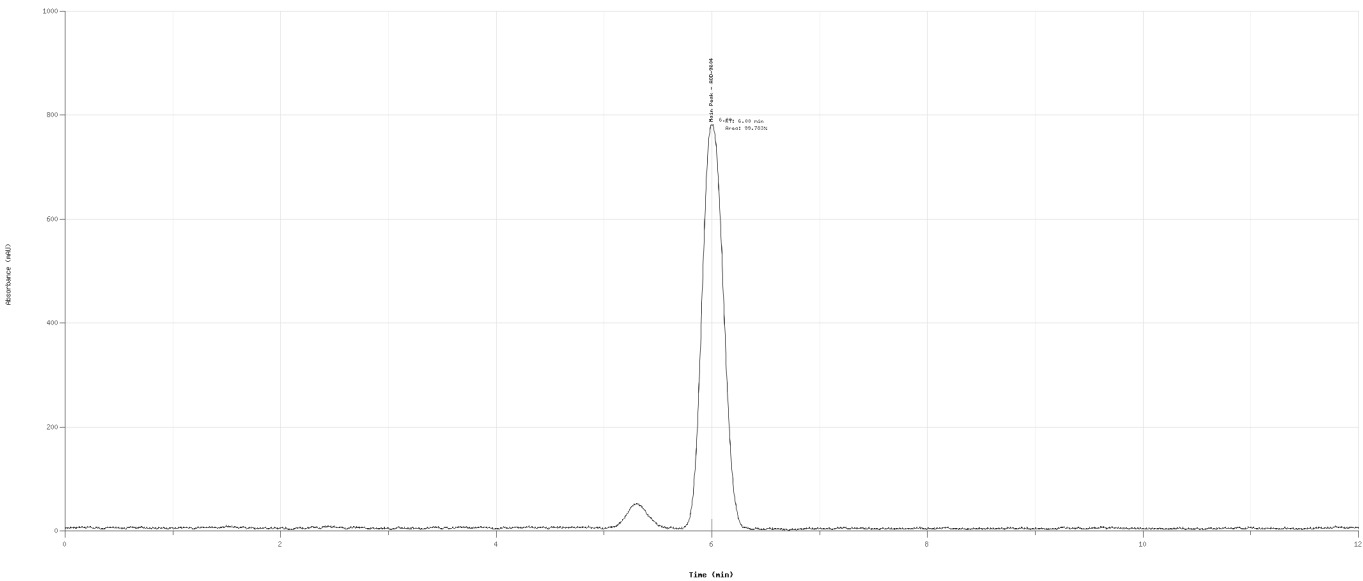

METHOD

TIME	H2O + 0.1% TFA	ACN + 0.1% TFA
0min	92%	8%
2min	92%	8%
10min	10%	90%
12min	10%	90%


TECHNICAL NOTE

Analytical conclusions are limited to the tested sample and the conditions described in this report. Total runtime: 12.0 minutes.

COMMENTS

The recorded analytical profile confirms that the sample meets the defined specifications for identity, purity and impurity profile using the applied RP-HPLC-UV method.

LC-MS MASS SPECTRUM

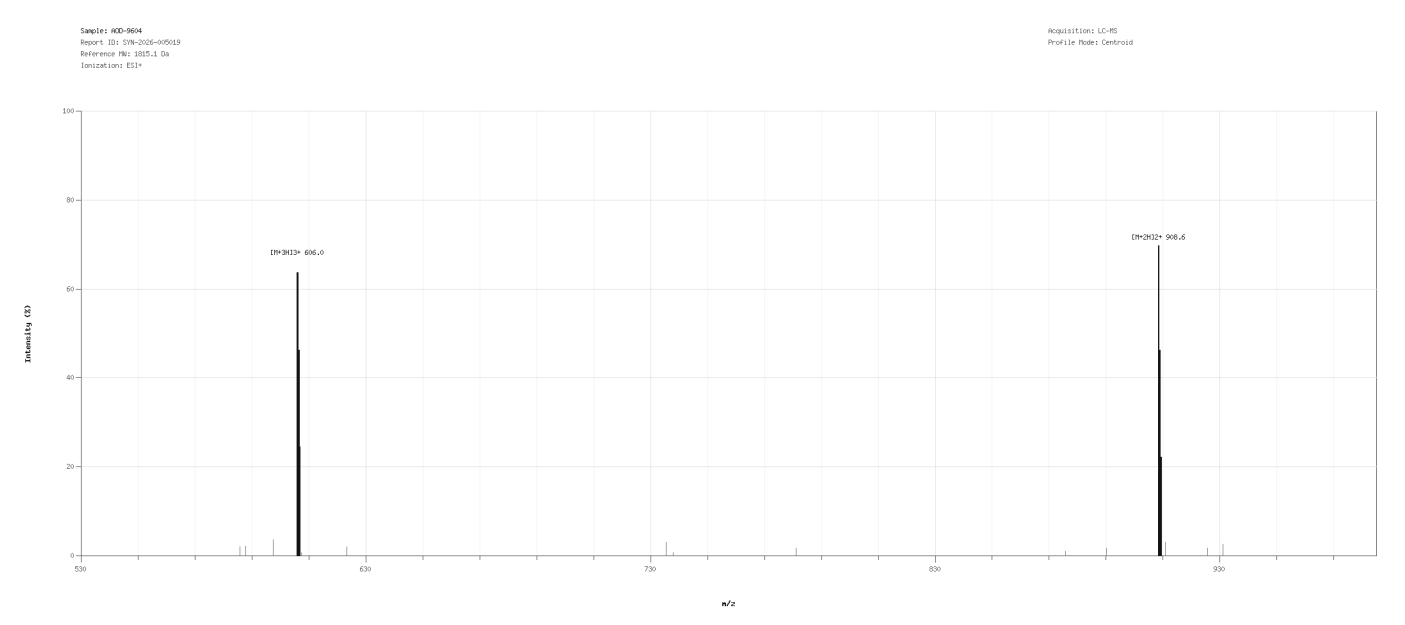

Charge-state distribution and isotope clustering are consistent with the expected analytical mass response.

ANALYSIS & METHODOLOGY

STANDARD PEPTIDE PURITY, MASS & ADDITIONAL TESTING

Sample analyzed by reverse-phase HPLC with UV detection using controlled analytical parameters. Chromatographic behavior and purity response were assessed through internal laboratory criteria.

BIOBURDEN

TEST	RESULT	UNIT	REPORTING LIMIT
Total Aerobic Microbial Count USP <61>/Eur. Ph. 2.6.12. Plate Count Method	Not detected	CFU/g	>= 1000
Total Yeast and Mold Count USP <61>/Eur. Ph. 2.6.12. Plate Count Method	Not detected	CFU/g	>= 100

ENDOTOXIN ANALYSIS

TEST	RESULT	UNIT	REPORTING LIMIT
Bacterial Endotoxin USP<85>/ Eur. Ph. 2.6.14. Bacterial Endotoxin Chromogenic Test	< 0.005	EU/mg	NMT 0.5

HEAVY METALS

TEST	RESULT	UNIT	REPORTING LIMIT
Arsenic Elemental Impurities Screening	Not detected	ppm	>= 1.5
Cadmium Elemental Impurities Screening	Not detected	ppm	>= 0.5
Cobalt Elemental Impurities Screening	Not detected	ppm	>= 25
Lead Elemental Impurities Screening	Not detected	ppm	>= 1.5

TEST	RESULT	UNIT	REPORTING LIMIT
Nickel Elemental Impurities Screening	Not detected	ppm	>= 25
Quicksilver Elemental Impurities Screening	Not detected	ppm	>= 1.5
Vanadium Elemental Impurities Screening	Not detected	ppm	>= 25

TECHNICAL APPENDIX

This appendix documents the analytical methodology, instrumentation and acceptance criteria applied for the evaluation of the sample.

COMPOUND REFERENCE

PARAMETER	AOD-9604
PUBCHEM CID	71300630
MOLECULAR FORMULA	C78H123N23O23S2
MOLECULAR WEIGHT	1815.1 g/mol

METHOD SPECIFICATION

PARAMETER	RP-HPLC-UV METHOD FOR PEPTIDE ANALYSIS
ANALYTICAL MODE	Purity assessment of peptide sample by RP-HPLC-UV
COLUMN	C18 peptide column, 150 x 4.6 mm equivalent
MOBILE PHASE A	Water + 0.1% TFA
MOBILE PHASE B	Acetonitrile + 0.1% TFA
FLOW RATE	0.8 mL/min
DETECTION	UV 220 nm
INJECTION VOLUME	10 uL
RUNTIME	12.0 min
SAMPLE DILUENT	Aqueous organic diluent compatible with peptide analysis
SAMPLE PREPARATION	Diluted, mixed and clarified before injection

INSTRUMENT PLATFORM

PARAMETER	STANDARD HPLC-UV PLATFORM
SYSTEM TYPE	Analytical HPLC system
DETECTOR	UV/VIS detector
ACQUISITION	Chromatographic acquisition and integration software
REVIEW MODE	Retention-time and response-profile review
WORKFLOW NOTE	Used for routine peptide purity and identity screening

ANALYTICAL CRITERIA

PARAMETER	ACCEPTANCE FRAMEWORK	BASIS
IDENTITY	Retention-time and profile agreement with reference expectations	Chromatographic identity review
PURITY	NLT 98.0% unless a stricter report-specific specification is declared	Integrated RP-HPLC-UV purity profile
QUANTITY	Measured content reviewed against the declared sample strength	Report-level analytical summary
BIOBURDEN	Not detected or within stated reporting limits	Microbial screening table
ENDOTOXIN	Below stated reporting limit / internal screening threshold	Endotoxin analysis table
HEAVY METALS	Below individual reporting limits where screened	Elemental impurities screening table

This appendix documents the analytical methodology, instrumentation and acceptance criteria applied for the evaluation of the sample.

VERIFICATION

Verify this report through the official Synaptica Analytics verification page using the details below.
Verification URL synaptica-labs.com/verify-report
Report ID SYN-2026-005019
Verification Key VK-PETQ-6VSP

SCAN TO VERIFY



DIGITAL SIGNATURE

A stylized, handwritten signature in black ink, appearing to read 'M. Saar'.

DIGITALLY SIGNED BY:
Martin Saar
Date: 2026.07.17
14:46:59 +02'00'
Director

info@synaptica-labs.com
Analysis date: Jul 15, 2026
Report generated: Jul 17, 2026
Analytical testing performed by Synaptica Analytics -
Analytical Services Division

Synaptica Analytics
SYN-2026-005019
Laboratory Analysis Report
VK-PETQ-6VSP

VERIFY AT
synaptica-labs.com/verify-report